

Self Food Quality Care

Consumption of vinegar with free mineral acid (fake vinegar) is health hazard. Vinegar must contain no free mineral acid. Testing of vinegar sample in laboratory takes time. The Department of Medical Sciences has developed kit for consumers to test free mineral acid, with fast and reliable result.

Health Impact

Consuming free mineral acid can cause severe abdominal pain, and gastritis.

Target Sample

Vinegar, vinegar with chili

Number of Test / Kit

100 Tests

Sensitivity

Detection limit 0.2%


Test Kit Tool

- 1 Bottle of free mineral acid testing powder
- 1 Plastic cup
- 1 Plastic spoon
- 1 Sheet of testing kit handbook

Procedure




1. Pour 5 ml. of vinegar into plastic cup

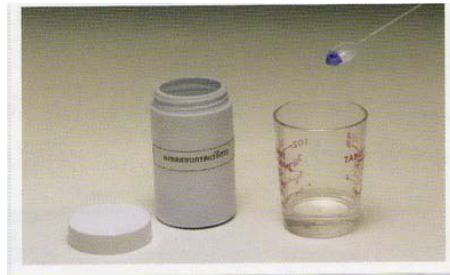

2. Add small amount of free mineral acid powder into plastic cup and shake until dissolve.

Evaluation

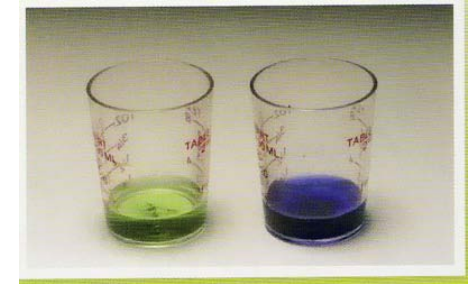

- If liquid in cup is purple, such sample is real vinegar (Consumable)
- If liquid in cup is blue or green, such sample is fake vinegar (unable to consume)

Test Kit for Free Mineral Acid in Vinegar (Fake Vinegar)



Place to Order

Asianmedic Co., Ltd

E-Mail : foodtest@asianmedic.com
www.asianmedic.com

Precaution

- Vinegar sample is acidity, when contacts skin, wash with water and soap.
- Keep test kit out of children reach.

Keeping / Age

Keep at room temperature / 2 years